

scatterBrain (V 1.230)

This document is a work in progress. Some scatterBrain features are yet to be documented, or have only been documented in a cursory manner.

Contents

Introduction	2
Installing scatterBrain	2
Preferred Installation	2
Alternative Installation	2
Upgrading scatterBrain	3
Starting scatterBrain	3
First Steps with scatterBrain	5
The Experiment File	5
The Live Log File	7
Loading Data	8
Plot Window.....	10
Profile tree	10
Subtracting Patterns	10
Axis scaling	11
Zoom in/out	11
Profile tree context menu	11
Profile Window Menu Functions	12
Save Profiles.....	12
Intensity Normalisation and Calibration	13
Q Calibration	18
Image Window	20
Main Program Menu.....	20
File.....	20
Acquire	20
Tools.....	20
Settings.....	20
Control Tabs.....	20
File Tree.....	20

Introduction

scatterBrain is a software package developed at the Australian Synchrotron SAXS/WAXS beamline for data acquisition, data reduction, and initial data analysis. It serves as the main user interface for scientists conducting experiments at the Australian Synchrotron SAXS/WAXS beamline. A take home version is available for users to complete their data analyses after their beamtime has finished.

scatterBrain is based on the SAXS15ID software originally written at the ChemMatCARS beamline of the Advance Photon Source Chicago. It is cross-platform compatible (Windows, Linux, OSX), and features a modern, drag-and-drop interface.

Installing scatterBrain

Preferred Installation

The preferred installation of scatterBrain is straight forward:

1. Download the Installation Package from:
<http://www.synchrotron.org.au/index.php/aussyncbeamlines/saxswaxs/software-saxswaxs>
This package is approximately 80 MB and contains a runtime distribution of IDL, compiled scatterBrain code (scatterBrainAnalysis.sav), a launcher, and various resource files.
2. Unarchive the package to any preferred location
3. Launch scatterBrain using the following launcher in the scatterBrainAnalysis directory:
 - a. scatterbrainanalysis.exe for Windows
 - b. scatterbrainanalysis .sh for Linux
 - c. scatterbrainanalysis .app for OS X
4. The first time scatterBrain is launched a hidden directory is created in your home folder. E.g., for windows this might be:
c:\documents and settings\stephen\.idl\australiansynchrotron\scatterbrain-###-windows\
The scatterBrainSettings.xml file is stored in this directory, and contains “user preference” settings, and a recent file list. This file persists even with upgrades to scatterBrain.

Alternative Installation

In rare cases scatterBrain doesn't launch when installed using the preferred method. Furthermore the “movie export” doesn't work using the preferred installation method due to licensing restrictions. The following install is recommend if you require movie mode or have other installation issues.

1. Complete steps 1 and 2 of the preferred installation method above
2. Create an account at www.exelisvis.com . Account creation is free.


- Download and install IDL from <http://www.exelisvis.com/Downloads/Download.aspx?product=IDL> . The download size is between 400 and 600 MB depending on your operating system. IDL can be downloaded for free. Precompiled code can be run without a license. A license is required only if you are writing your own IDL code.
- Launch scatterBrainAnalysis using scatterBrainAnalysis.sav – note that this is a different file to the one used in the preferred installation method.
- The scatterBrainSettings file mentioned above will also be created with this installation method.

Upgrading scatterBrain

Once installed, scatterBrain can be upgrade through the scatterBrain interface itself. In future versions you will be alerted of updates at started up. However, currently you manual check for updates using the option under the help menu. The update procedure usually just replaces the scatterBrainAnalysis .sav file which contains the actual code, however, larger updates may require the resource files to be updated or a new IDL distribution downloaded. You will be alerted of these larger updates through the scatterBrain interface and/or website.

Starting scatterBrain

Upon launching scattering you will be greeted with a window not unlike this (you won't see this window if you have used the alternative installation method):

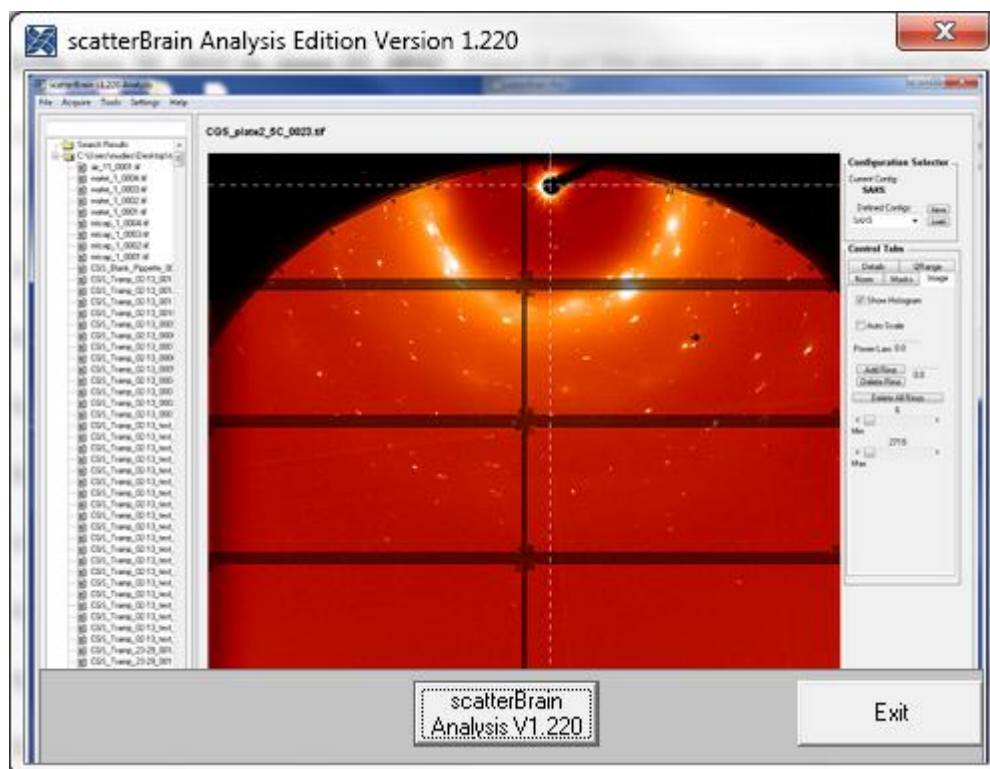


Click on the scatterBrain Analysis button at the bottom centre to continue. You will then see yet another splash screen:




Click anywhere on this splash screen to continue. This second splash screen appears because we are using the IDL virtual machine. The only way to get rid of this is to buy an IDL license.

scatterBrain will then open, as two windows. The image window:

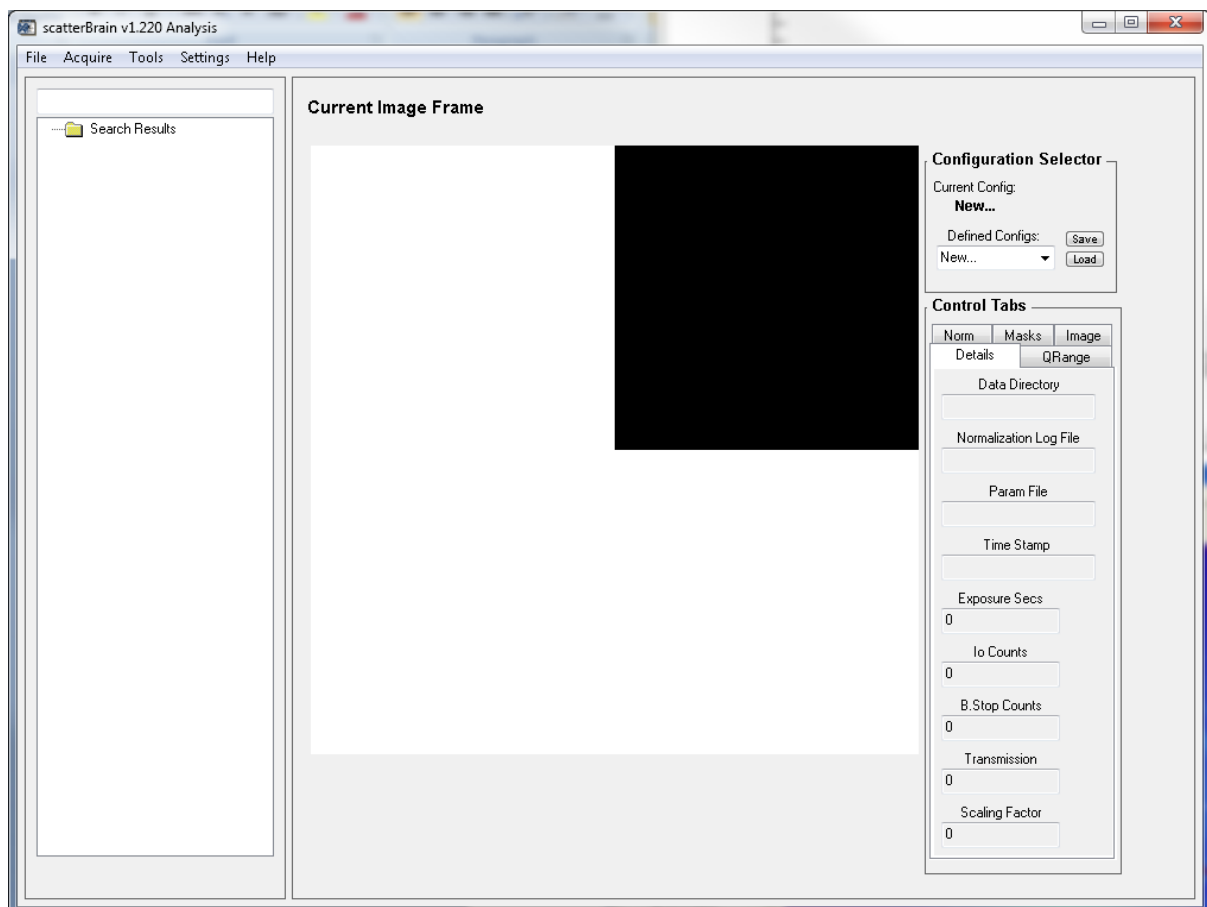

and the plot window:

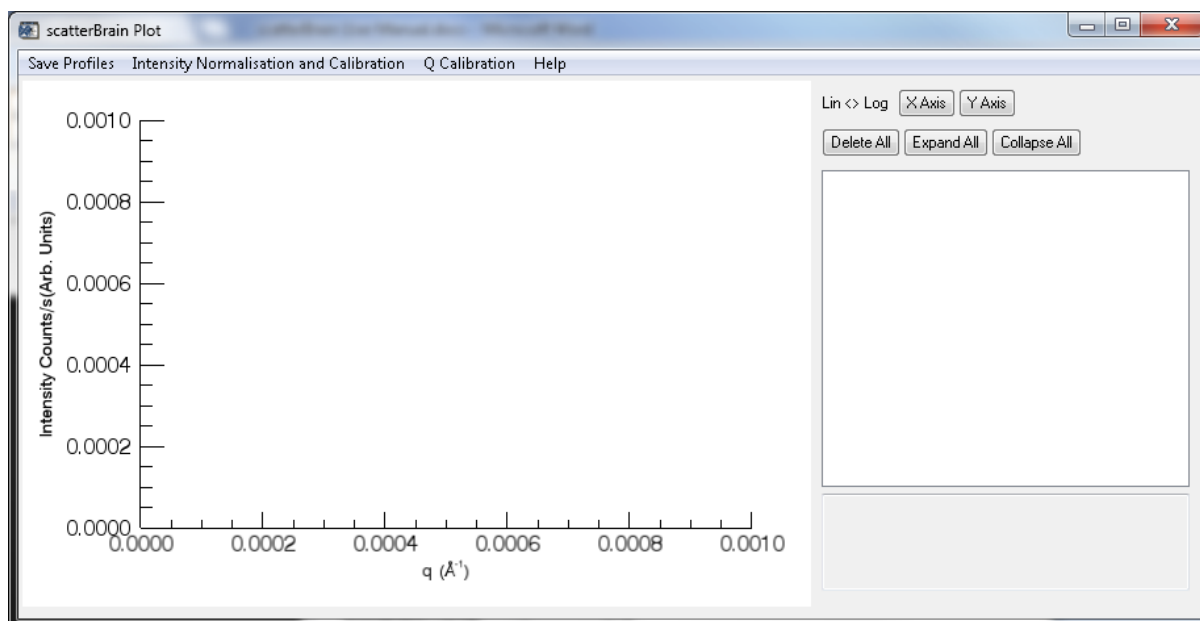


All screenshots in this document are from the Windows 7 version of scatterBrain, because it looks prettiest on that platform. scatterBrain looks terrible on Linux/OS X, but that is underlying issue with IDL that I hope they upgrade soon. Sigh...

First Steps with scatterBrain

The Experiment File

To start looking at data with scatterBrain, you first must open a scatterBrain experiment. To load the experiment file:

File->Get Experiment

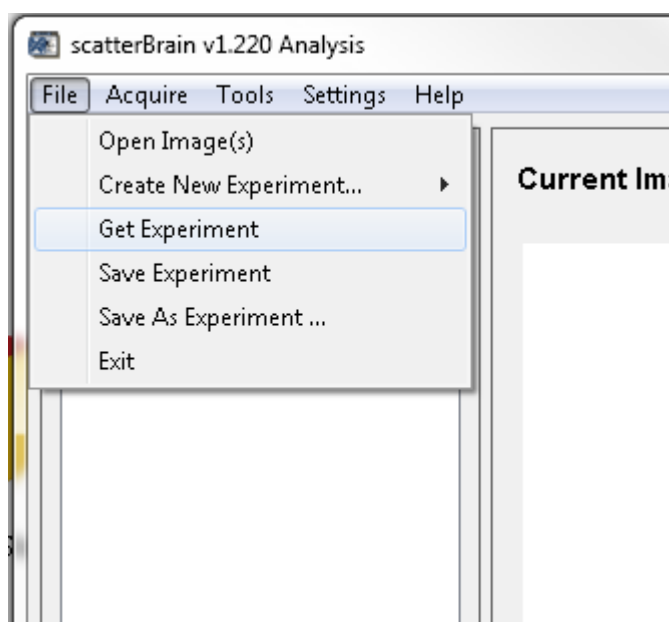

An experiment file is in the XML format, and looks something like this:

```

1  > Hawley_4329_SAXS_day2.xml
2  <!-->
3  <!-->
4  <!-->
5  <!-->
6  <!-->
7  <!-->
8  <!-->
9  <!-->
10 <!-->
11 <!-->
12 <!-->
13 <!-->
14 <!-->
15 <!-->
16 <!-->
17 <!-->
18 <!-->
19 <!-->
20 <!-->
21 <!-->
22 <!-->
23 <!-->
24 <!-->
25 <!-->
26 <!-->
27 <!-->
28 <!-->
29 <!-->
30 <!-->
31 <!-->
32 <!-->
33 <!-->
34 <!-->
35 <!-->
36 <!-->
37 <!-->
38 <!-->
39 <!-->
40 <!-->
41 <!-->
42 <!-->
43 <!-->
44 <!-->
45 <!-->
46 <!-->
47 <!-->
48 <!-->
49 <!-->
50 <!-->
51 <!-->
52 <!-->
53 <!-->
54 <!-->
55 <!-->
56 <!-->
57 <!-->
58 <!-->
59 <!-->
60 <!-->
61 <!-->
62 <!-->
63 <!-->
64 <!-->
65 <!-->
66 <!-->
67 <!-->
68 <!-->
69 <!-->
70 <!-->
71 <!-->
72 <!-->
73 <!-->
74 <!-->
75 <!-->
76 <!-->
77 <!-->
78 <!-->
79 <!-->
80 <!-->
81 <!-->
82 <!-->
83 <!-->
84 <!-->
85 <!-->
86 <!-->
87 <!-->
88 <!-->
89 <!-->
90 <!-->
91 <!-->
92 <!-->
93 <!-->
94 <!-->
95 <!-->
96 <!-->
97 <!-->
98 <!-->
99 <!-->
100 <!-->
101 <!-->
102 <!-->
103 <!-->
104 <!-->
105 <!-->
106 <!-->
107 <!-->
108 <!-->

```

Zooming in on the important bits:

```

79 <scatterBrain>
80
81 <Parameters>
82 <PVMap>
83 <PV ACQUIRE="false" ACQUIREHIGH="10" ACQUIRELOW="-10" LOG="false" READ="SR13ID01HU02IOC01:SMPL_TBL_X_MTR.RBV" SET="SR13ID01HU02"
84 </PVMap>
85 <DetectorMap>
86 <DETECTORDEF AUTOLOAD="1" BASEPV="13PIL1:" CAMPV="cam1:" CONTROL="1" FILEPV="" IMAGEPV="image1:" PIXELSIZE="
87 <DETECTORDEF AUTOLOAD="0" BASEPV="13PIL2:" CAMPV="cam1:" CONTROL="1" FILEPV="" IMAGEPV="image1:" PIXELSIZE="
88 </DetectorMap>
89 <MaskMap>
90 <USERMASK AUTO="false" MASKNAME="Pilatus 1M Frame" SHAPE="Polygon">[487.000,487.000,0.000000,0.000000,487.000,487.000,0.000000,0
91 <USERMASK AUTO="false" MASKNAME="Mask 2" SHAPE="Polygon">[-2.00000,28.0000,123.000,247.000,391.000,788.000,837.000,929.000,985.
92 <USERMASK AUTO="false" MASKNAME="Mask 3" SHAPE="Polygon">[-2.00000,56.0000,-3.00000][966.000,1044.00,1045.00]</USERMASK>
93 <USERMASK AUTO="false" MASKNAME="Mask 4" SHAPE="Polygon">[838.000,843.000,843.000,839.000][862.000,862.000,858.000,858.000]</US
94 <USERMASK AUTO="false" MASKNAME="Mask 5" SHAPE="Polygon">[674.000,680.000,680.000,675.000][86.0000,86.0000,82.0000,81.0000]</US
95 </MaskMap>
96 <CONFIGURATION NAME="SAXS">
97 <CAMERADEFs BEAMX="522.000" BEAMY="25.0000" DETECTOR="Pilatus 1M" LENGTH="1028.11" STOPACTIVE="1" STOPAI
98 <USERMASKS>[Pilatus 1M Frame,Mask 2,Mask 3,Mask 4,Mask 5][1,1,1,1,1]</USERMASKS>
99 <COUNTERDEFs BEAMSTOP="IBS" INCIDENT="I0" TRANSMISSION="IT"></COUNTERDEFs>
100 <NORMALISATION ABSICAL="0.000000" I0NORM="0.000000" IBSNORM="0.000000" NORMTYPE="0" USEABSICAL="0"></I
101 </CONFIGURATION>
102 </Parameters>
103
104 <Experiment>
105 </Experiment>
106
107 </scatterBrain>
108

```

```

IOC01:SMPL_TBL_X_MTR.RBV" TRANS="false" TRANSHIGH="-8" TRANSLOW="-9"></PV>

0.17200000" SOFTWARETRIGGER=" 1" XSIZE=" 981" YSIZE=" 1043">Pilatus 1M</DETECTORDEF>
0.17200000" SOFTWARETRIGGER=" 0" XSIZE=" 981" YSIZE=" 195">Pilatus 200K</DETECTORDEF>

0.000000,487.000,487.000,0.000000,0.000000,487.000,487.000,0.000000,0.000000,487.000,487.000,494.000,494.000,981.000,981.000,494.000,4
.000,983.000,-4.000000][286.000,225.000,123.000,33.0000,-2.00000,6.00000,23.0000,82.0000,132.000,-2.00000,-3.00000]</USERMASK>

SERMASK>
SERMASK>

ANGLE=" -25.0000" STOPOFFSETX=" 0.000000" STOPOFFSETY=" 0.000000" STOPRADIUS=" 15.0000" STOPWIDTH=" 17.0000" WAV

/NORMALISATION>

```

This file contains all the configuration parameters required to generate the normalised 1D profiles, such as beamcentre, camera length and wavelength, masking, absolute intensity calibration, pixel size, etc. It is typically created at the start of the experiment by the beamline scientists, although it can be created at a later date, with the parameters determined from calibration shots. It is overwritten whenever it is saved inside scatterBrain, there is no automatic saving. To save:

File->Save Experiment

You can also save the experiment to a new file using:


File->Save As Experiment

Initially the experiment file doesn't contain any information about the SAXS exposures (during an experiment this is stored in the live log, as explained in the next section), however, whenever the experiment is saved, the currently loaded log information is also saved into the experiment file. Hence, if you save the experiment with the full log file for your experiment loaded, in future you can forget about the log file. The experiment file also contains summed file normalisation parameters and detector parameters for running experiment.

The Live Log File

The Live log file is generated during data acquisition. It has an XML like format. Data is only ever appended to the file, it is never overwritten by the control system. A new line is added to the file with ever exposure, and many parameters are logged, but the most important is normalisation data.

The file looks like this:

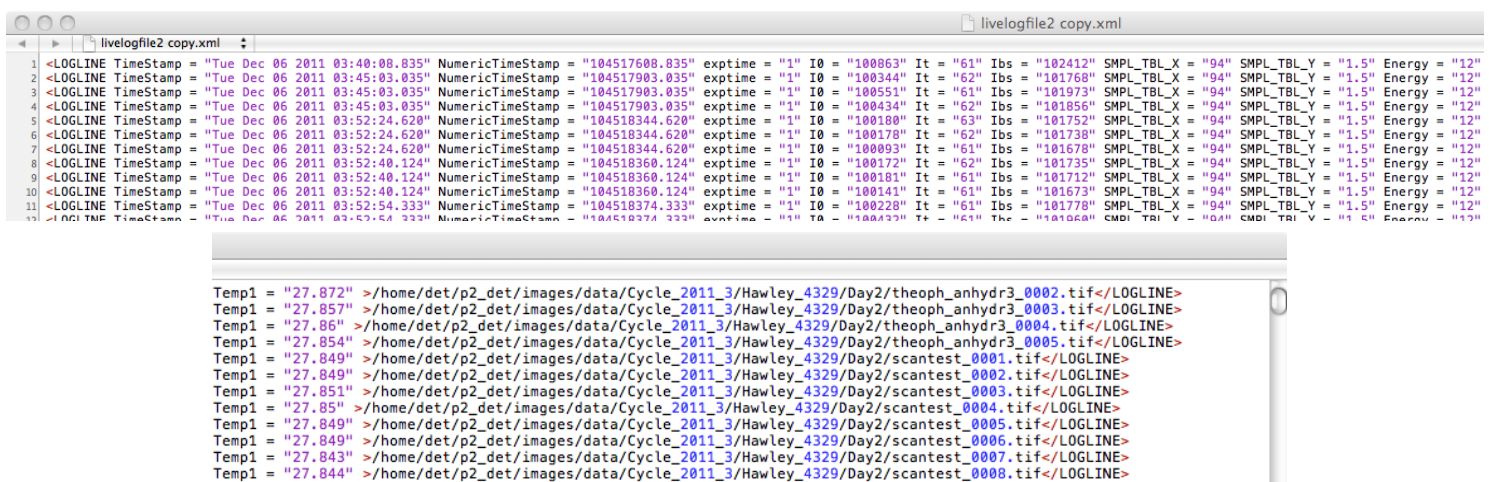

```


livelogfile2 copy.xml
livelogfile2 copy.xml
1 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:40:08.835" NumericTimeStamp = "104517608.835" exptime = "1" I0 = "100863" It = "61" Ibs = "102412" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
2 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:45:03.035" NumericTimeStamp = "104517903.035" exptime = "1" I0 = "100344" It = "62" Ibs = "101768" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
3 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:45:03.035" NumericTimeStamp = "104517903.035" exptime = "1" I0 = "100551" It = "61" Ibs = "101973" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
4 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:45:03.035" NumericTimeStamp = "104517903.035" exptime = "1" I0 = "100434" It = "62" Ibs = "101856" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
5 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:24.620" NumericTimeStamp = "104518344.620" exptime = "1" I0 = "100180" It = "63" Ibs = "101752" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
6 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:24.620" NumericTimeStamp = "104518344.620" exptime = "1" I0 = "100178" It = "62" Ibs = "101738" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
7 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:24.620" NumericTimeStamp = "104518344.620" exptime = "1" I0 = "100093" It = "61" Ibs = "101678" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
8 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:40.124" NumericTimeStamp = "104518360.124" exptime = "1" I0 = "100172" It = "62" Ibs = "101735" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
9 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:40.124" NumericTimeStamp = "104518360.124" exptime = "1" I0 = "100181" It = "61" Ibs = "101712" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
10 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:40.124" NumericTimeStamp = "104518360.124" exptime = "1" I0 = "100141" It = "61" Ibs = "101673" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
11 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:54.333" NumericTimeStamp = "104518374.333" exptime = "1" I0 = "100228" It = "61" Ibs = "101778" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"
12 <LOGLINE TimeStamp = "Tue Dec 06 2011 03:52:54.333" NumericTimeStamp = "104518374.333" exptime = "1" I0 = "100423" It = "61" Ibs = "101868" SMPL_TBL_X = "94" SMPL_TBL_Y = "1.5" Energy = "12"

Temp1 = "27.872" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/theoph_anhydr3_0002.tif</LOGLINE>
Temp1 = "27.857" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/theoph_anhydr3_0003.tif</LOGLINE>
Temp1 = "27.86" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/theoph_anhydr3_0004.tif</LOGLINE>
Temp1 = "27.854" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/theoph_anhydr3_0005.tif</LOGLINE>
Temp1 = "27.849" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0001.tif</LOGLINE>
Temp1 = "27.849" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0002.tif</LOGLINE>
Temp1 = "27.851" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0003.tif</LOGLINE>
Temp1 = "27.85" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0004.tif</LOGLINE>
Temp1 = "27.849" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0005.tif</LOGLINE>
Temp1 = "27.849" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0006.tif</LOGLINE>
Temp1 = "27.843" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0007.tif</LOGLINE>
Temp1 = "27.844" >/home/det/p2_det/images/data/Cycle_2011_3/Hawley_4329/Day2/scantest_0008.tif</LOGLINE>

```

As you can see each line has a timestamp, exposure time, normalisation parameters, information about sample table position, etc, along with the filename of the image collected. This file doesn't contain information about the camera setup, etc, as this is also stored in the experiment file as described above. To load the live logfile:

Acquire->Select Live Log

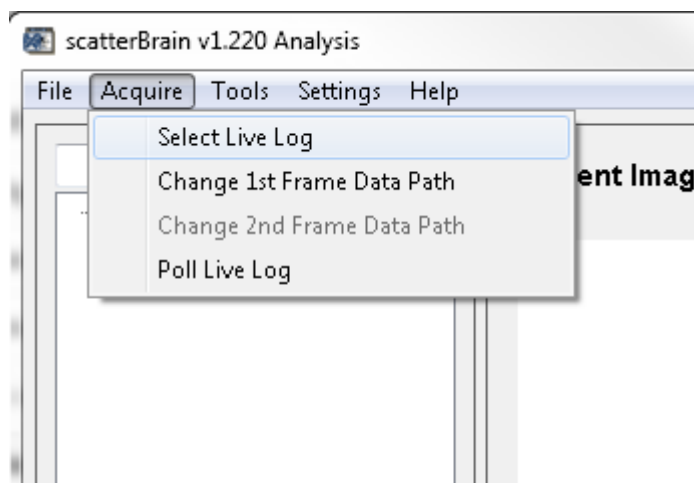

During data acquisition at the beamline this file is being constantly updated, and scatterBrain will display these updates in real time. Anytime you save your experiment, the current state of the live log will also be saved into your experiment file. This means if you save your experiment file once you have collected all of your data, or if at a later date you open your experiment file and log file, and then save your experiment file, that the experiment file will then contain all of the information about your experiment. From that point forward you will not need to open the live log when you open your experiment.

Loading Data


After opening your experiment file (and possibly your live log file) the file list tree on the left hand side of the image window will be populated with a list of all the files known by the experiment (and/or log file). To load an image you:

double-click it


or

drag and drop it into the file tree on the Plot window

The 2D image will be displayed on the main screen and the 1D azimuthally averaged plot will display on the plot screen. The image window will now look something like the following. Note the file list on the left:

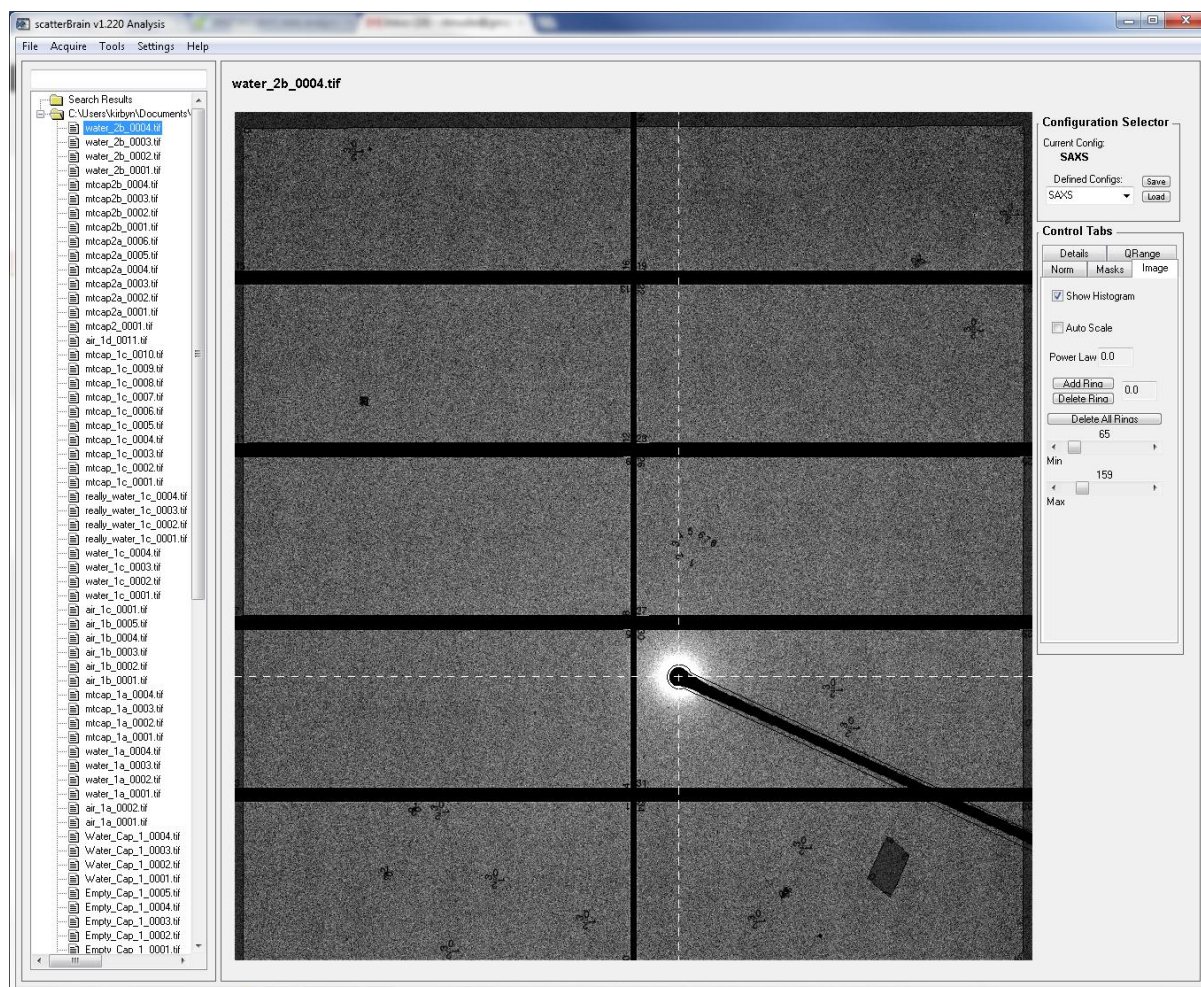

and after loading two images, the plot window may look like this:

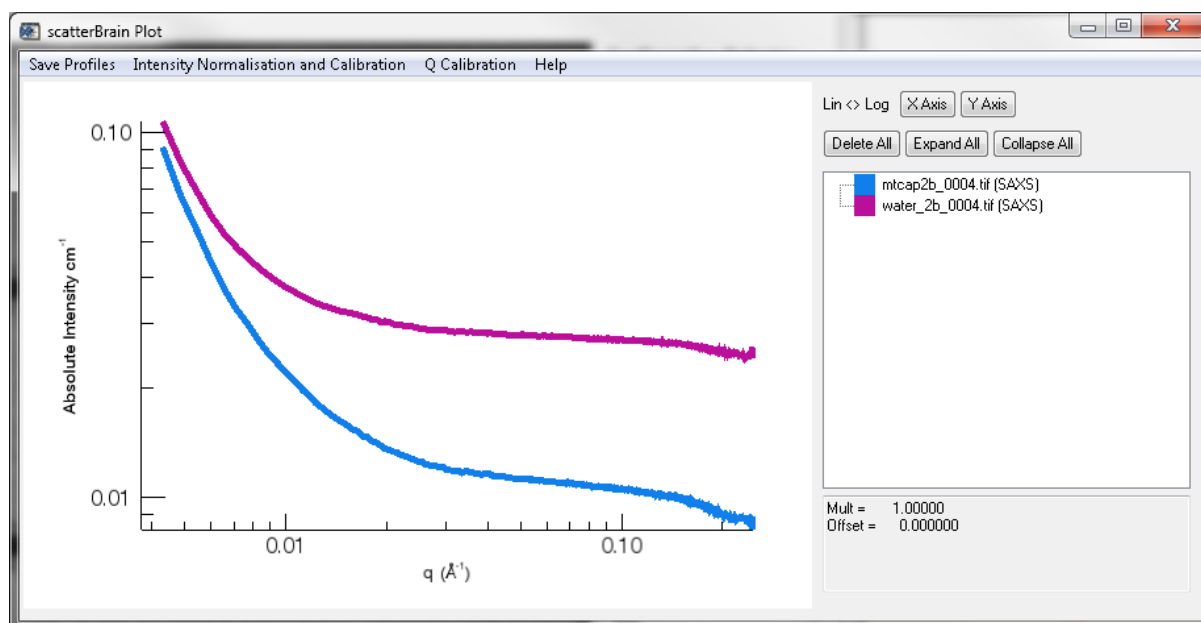

Plot Window


The plot window has many features for displaying, using and generating 1-dimensional patterns, and is the window where most 1-d pattern operations are performed.

Profile tree

The profile tree is a list of all the 1-D patterns currently held in memory, and is shown on the right hand side of the plot screen.

Subtracting Patterns

Profile subtraction is done by drag and drop: DROP SAMPLE PROFILES ONTO BACKGROUND PROFILES. You can do multiple patterns at once by selecting them with shift-left click (for continuous files in the tree) or CTRL-left click (for discontinuous groups of files in the tree). The background

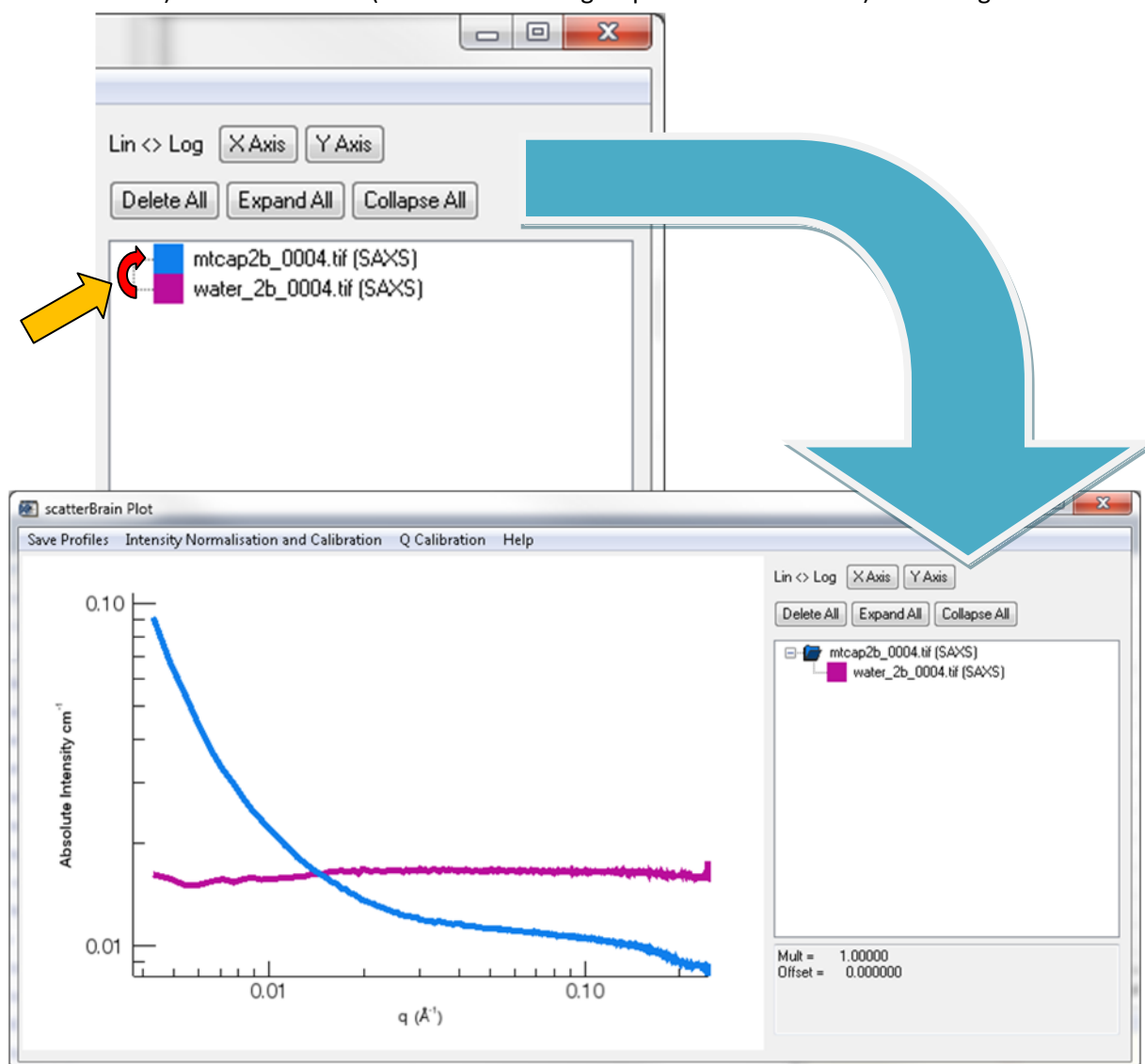

pattern will stay at top level hierarchy in the profile tree and the subtracted sample pattern(s) will move to second level in the tree hierarchy (i.e. the sample profiles moving slightly to the right). The background will also now show a “-” or “+” symbol which allows the sample profiles under each background to be shown or hidden in the tree listing. The plot window will now display the subtracted

profile(s). Background profiles will still be displayed by default, although you can hide individual profiles (either samples or backgrounds) using the profile tree context menu.

To unsubtract a profile in the profile tree, just drag it off its background profile to anywhere in the blank (white) area of the profile tree. To change which background is subtracted from a sample profile, drag the sample profile from one background to another in the profile tree.

Plot highlighting: If you click on one (or more) plots in the tree display, selected plots(s) will be highlighted with a thicker line style. This is a fast way to see which pattern is which on the main plot window when there are many pattern overlaid on screen. You can also move up/down in the plot tree using a keyboard (up/down keys) which can be more convenient than using a mouse.

Axis scaling

Linear/log axes. You can toggle the X and Y axes between linear and logarithmic scales using the buttons above the plot tree.

Zoom in/out

Left-drag a box on the main plot area, and the plot will zoom in to that area of the plot. To un-zoom, right click in the plot area and the plot will unzoom and autoscale. When new plots are added to the plot window, the axes hold their current scale ranges and do not autoscale: if the new plot doesn't fit the current scale (e.g. you can't see it at all) then right-click to autoscale.

Profile tree context menu

The context menu of the plot tree contains many functions which can act on existing profiles. Select the profile(s) you want to use with left-click, shift-left click or ctrl-left click, then

Right-Click anywhere in the plot tree

and the context menu will appear. Most actions are self explanatory, so just some brief guidance:

- **Delete** removes the profile from the profile tree, the plot area and from memory. It does not delete the original data (tif file). The "**Delete All**" button above the profile tree will delete all profiles.
- **Apply Parameters** recomputes the profile, in order to apply any changes in camera parameters
- **Multiply** and **Offset** allow you to apply multiplication or additive correction to a profile. The defaults are 1 and 0 respectively (i.e. no correction). If you have small subtraction errors these can be useful. Each correction is applied to whichever pattern(s) is selected before you apply the correction – normally you would want to only correct one profile at time. Any multiplication or addition terms for the currently selected profile are displayed below the profile tree.
- **Fit Peak** fits a linear combination of a Gaussian peak and a linear background to the part of the currently selected profile (selected in the profile tree) which is displayed in the plot area. To fit a peak well you should normally zoom in on the peak first. Basic properties of the peak fit (position, width etc.) are temporally displayed at the bottom of the plot area.
- **Hide, Show** and **Show<>Hide** allow you to hide a profile from the plot area without removing it from the profile tree, to re-display a hidden profile, or to toggle from one state to another.
- **Change Colour** allows the profiles colour to be changed. A colour table dialog box appears for handling all plot colours. This is explained latter in this document.

- **Line Width** and **Opacity** alter the line properties on the plot area.
- **Protein Tools** allows some basic protein-specific data analysis to be run. Currently this allows **Autoprod** (ATSAS suite) to be run on one profile at a time (the last one selected). The fit data is displayed on screen and a temporary data file is written to the operating system's temporary directory.

Profile Window Menu Functions

Save Profiles

This menu item has two options for saving profiles to disk in a space delimited ASCII format.

The first option

Save ALL profiles to ONE large file

saves all the currently plotted profiles into one single file. There are three columns of data for each profile, q vector, Intensity, Error in Intensity. The filenames for the each of the profiles is written above each set of data. When this menu item is selected a dialog box appears in order to select the file into which the data will be saved. The second option

Save EACH file to individual profiles

Saves a separate file for each profile plotted. These files have three columns each, q vector, Intensity Error in Intensity. When you select this option a dialog box appears for selecting the directory for output. The profiles will be saved into the files with the same name as the original data, with a .dat file extension.


Note that both of these export functions honour any multiplicative and offset factors applied, as well as background subtraction.

Intensity Normalisation and Calibration


1. Delete all data from the plot window.
2. From q-calibration, check the camera length is set correctly. The typical camera lengths for proteins are:

- 959 mm
- 1576 mm
- 3349 mm
- 7200 mm.


Check photon energy is correct. We now use 11 000 eV: 12 keV is so last-year. If you need to change the camera length, type the correct length into both the "Calibrated Camera Length" and the "Current Camera Length" fields. If you make any changes, click "Apply Q Calibration", then close the window.

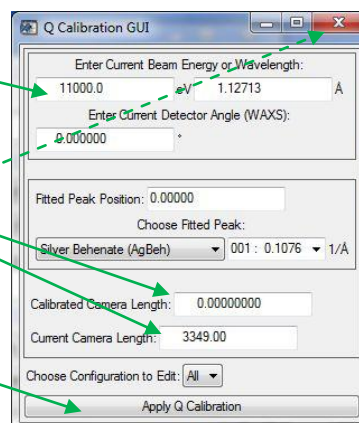

3. Load your airshot

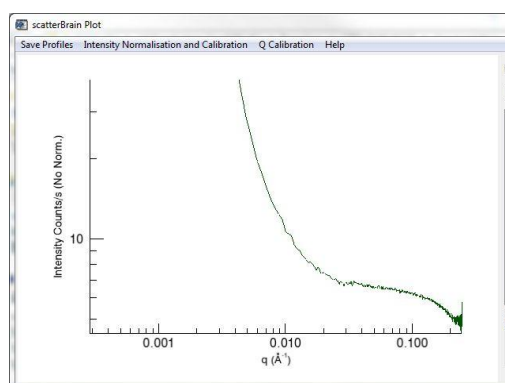


4. Open the Intensity Normalisation and Calibration menu.

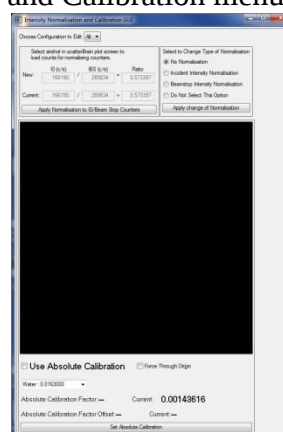


- Left click in the white space here


- Then select (click on) the air pattern here

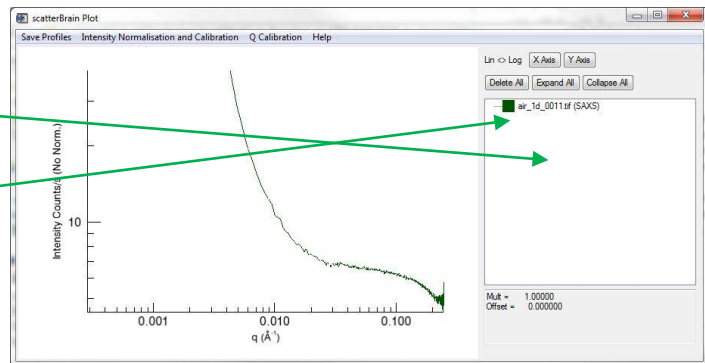

- This will update the "New" Io, IBS and Ratio normalisation inputs.

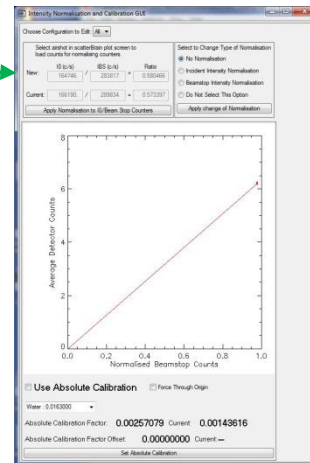

- Then click "Apply Normalisation to IO/Beam Stop Counters" to make the new numbers "Current". Both rows of numbers need to match.

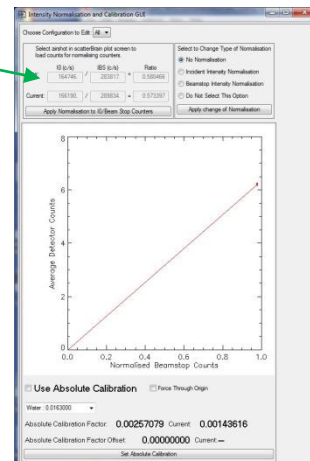


- Select the "Beamstop Intensity Normalisation" radio button


- Click Apply change of Normalisation button.

This will allow patterns to be normalised for incident and transmitted flux variation using the beamstop detector.


- Close the window.

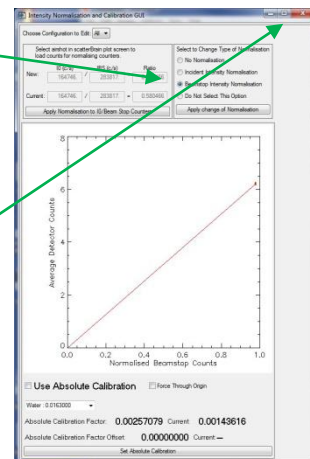

12. Delete the air pattern, and load the mtcap and water shots.

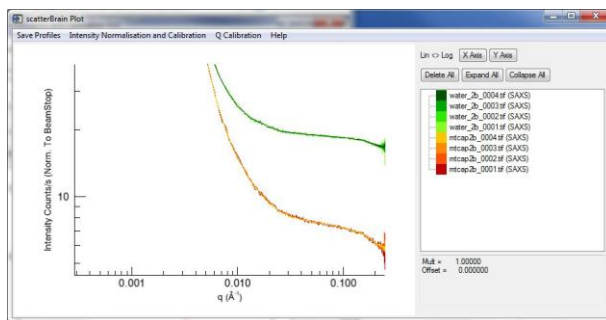


13. Subtract the mtcap shots from the water shots (pairing up 1,2,5, and 10 second exposures).

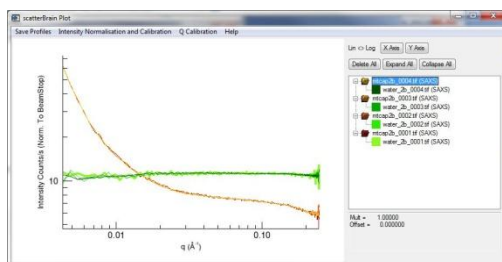

14. Load the "Intensity Normalisation and Calibration" Screen




15. Click in the white space

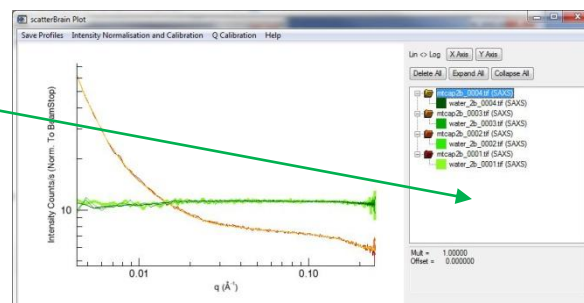

16. select all of the water shots (do not select the mtcap shots)


17. Zoom in (left-drag box on plot screen) to select flat part of water calibration (if needed).

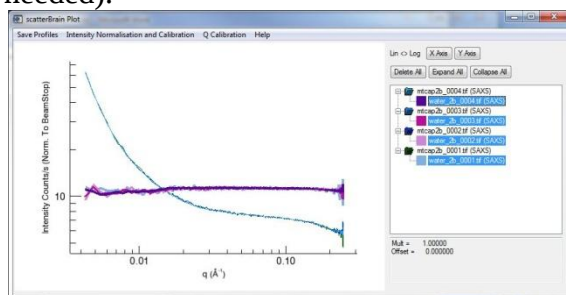

18. You should see a plot that looks like this.


19. Click on "Set Absolute Scaling".

The tick box next to "Use Absolute Calibration" will come on. A small scale factor offset (eg between ± 0.1 is OK). In V1.61, ignore status of normalisation mode at this point in the procedure (don't change it). There is a bug in the code making the current norm state not display correctly just at this part of the calibration sequence which will be eliminated in next version.

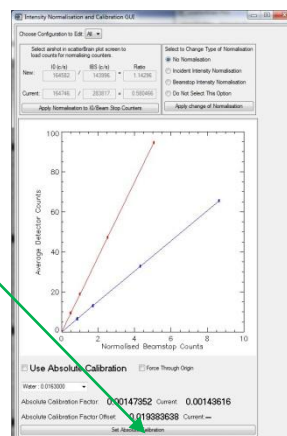


20. Close the window.

21. Delete all the patterns in the plot window and re-load air, mtcap and water shots. Subtract mtcap from water shots, and check you are happy with the results. The flat part of the water pattern (when its mtcap background is subtracted) must be close to 0.0163 cm^{-1} . If not, repeat software calibrations. If this is still wrong then you may need to re-collect data. You are also looking for an acceptably clean subtraction at low q, i.e. the subtracted water pattern is flat and doesn't tail up or down too much - how much depends on what scattering intensity you are then going to use, which depends on protein size, scattering power and concentration. The following is an good calibration data set at 3.3 m camera length:

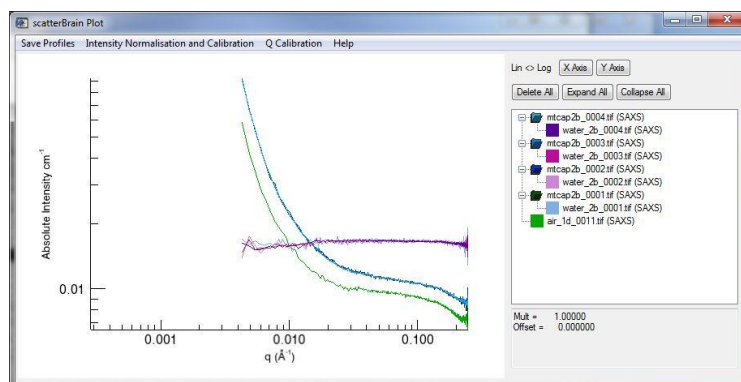

Q Calibration

This menu item opens a new window for energy/wavelength/camera length calibration.

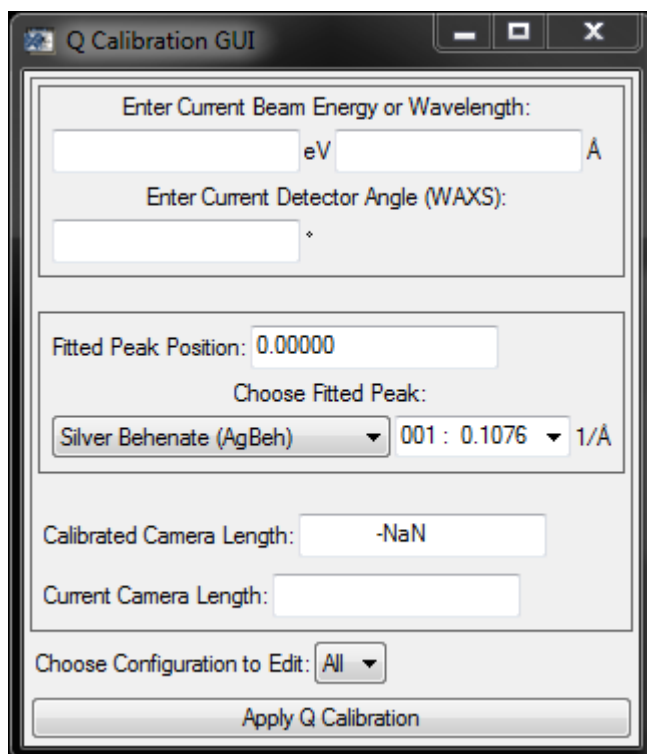

The beam energy or wavelength is entered/alterd with the keyboard (these fields are linked, so you only have to fill in one).

To use a calibration standard to set the camera length, you first load a measured pattern, zoom in to just a single peak in the plot area, right-click on that profile in the profile tree and select:

Fit Peak

The **Fitted Peak Position** will be updated in the Q Calibration GUI, and then the camera length will be calculated based on the reference standard (e.g. Silver Behenate, rat tail tendon, LaB6 or a user-configured standard) and the diffraction peak order currently selected from the drop-down lists. You can change standards and/or peaks as needed and the camera length will be re-calculated. The calculated (i.e. calibrated) camera length is shown in the "Calibrated Camera Length" box, and the actual camera length currently being used for q-scaling is shown in the "Current Camera Length". To apply the new calibration (to the Configuration(s) selected in the configuration drop-down (often just use "All") click the

Apply Q Calibration

button. This will update the camera length used for q-scaling and close the Q-calibration window. If you wish to use this window to only change the energy, without changing the camera length, make sure that the

Calibrated Camera Length

Is set to be the same as the

Current Camera Length

before applying.

Image Window

Is the nerve centre of scatterBrain. Not only does it display the image data, but includes the File Tree, tools for determining beamcentre, setting masks, as well as displaying the main program menu.

Main Program Menu

The main program menu has several

File

Acquire

Tools

Several tools are accessible from this menu item. The

Tools->Export Image->Export Current Image

Command exports the currently loaded image into a format that is convenient for viewing and printing. Typical file formats such as JPEG and PNG are offered. The exported file has the same colour table and scaling applied as the image in the image window. Note that these exported images have much less information encoded in them as the original data, as the bit depth is reduced from the native 32bits in order to view on screen or print.

Tools->Export Image->Export Current Image with Annotations

Is the same as above but also includes masking and beamcentre annotations.

The

Convert SAXS15ID log/saxs files to scatterBrain experiment

Allows data generated using SAXS15ID to be converted for use in scatterBrain.

Settings

Settings->General Settings

Command opens a window with many general settings available for scatterBrain. These include

- ***Pixel Bins*** – The number of pixel widths per histogram bin. Typically this should be set to 1 or 2.
- ***Number of Sectors*** – The number of sectors to use in a sector integration.

Control Tabs

The control tabs on the right hand side of the image window, enable access to various functions and information.

File Tree